WORLD INTELLECTUAL PROPERTY ORGANIZATION International Bureau



INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification 5: G01N 33/569, 33/68, A61K 39/21 C12P 21/08

(11) International Publication Number:

WO 92/15885

• | /

(43) International Publication Date:

17 September 1992 (17.09.92)

(21) International Application Number:

PCT/US92/01749

(22) International Filing Date:

6 March 1992 (06.03.92)

(30) Priority data:

668,266

11 March 1991 (11.03.91) US

(81) Designated States: AT (European patent), AU, BE (European patent), BR, CA, CH (European patent), DE (European patent), ES (European patent), FI, FR (European patent), GB (European patent), GR (European patent), IT (European patent), JP, KR, LU (European patent), MC (European patent), NL (European patent), NO, SE (European patent).

,

Published


With international search report.

(71) Applicant: IDEC PHARMACEUTICALS CORPORA-TION [US/US]; 11099 North Torrey Pines Road, Suite 160, La Jolla, CA 92037 (US).

(72) Inventor: KANG, Chang-Yuil; 1736 Buttercup Road, Encinitas, CA 92024 (US).

(74) Agents: DAHL, Nancy, K. et al.; Lyon & Lyon, 611 West Sixth Street, 34th Floor, Los Angeles, CA 90017 (US).

(54) Title: METHODS FOR SELECTING ANTIBODY REAGENTS; ANTI-IDIOTYPE ANTIBODIES; AND AIDS VACCINE FORMULATIONS

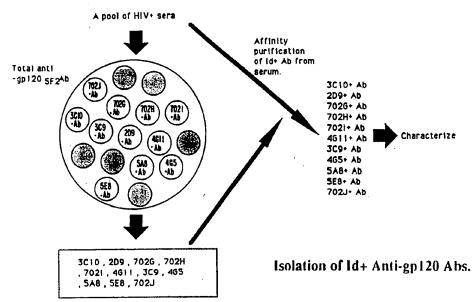

(57) Abstract

A panel of monoclonal anti-idiotype antibodies

A novel method is disclosed for selecting antibody reagents for vaccine formulations that are capable of inducing group-neutralizing anti-HIV antibody responses in humans. Also disclosed are anti-idiotype antibodies and hybridoma cell lines capable of producing such antibodies. Further, anti-idiotype antibody-based vaccines and formulations thereof are disclosed. The methods, antibodies, hybridomas and vaccine formulations of the present invention may be useful in the treatment and prevention of HIV infection.

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

America	¥1	binland	MI	Mali
• • • • • • • • • • • • • • • • • • • •	-		MN	Mongolia
			MR	Mauritania
Barbados				Malawi
Belgium	_	•		Netherlands
Burkina Faso	GN	Guinea		
Bulgaria	GR	Greece		Norway
•	HU	Hungary	PL	Poland
	IE	Ireland	RO	Romania
		Italy	RU	Russian Federation
		•	SD	Sudan
	_		SE	Sweden
	KP	•		Senegal
Switzerland		•	-	Soviet Union
Côte d'Ivoire			_	Chad
Cameroon	LI	Liechtenstein		
Czechostovakia	LK	Sri Lanka		Togo
	เม	Luxembourg	US	United States of America
•	MC	Monaco		
	Burkina Faso Bulgaria Benin Brazil Canada Central African Republic Congo Switzerland Côte d Ivoire Canicroon Czechislovakia Germany Denmark	Australia FR Barbados GA Belgium GB Burkina Faso GN Bulgaria GR Benin HU Brazil IE Canada IT Cantral African Republic JP Congo KP Switzerland Cöte d'Ivoire KR Canieroon LI Czechostovakia LK Germany LU Denmark MC	Australia FR France Barbados GA Gabon Belgium GB United Kingdom Borkina Faso GN Guinea Bulgaria GR Greece Benin HU Hungary Brazil IE Ireland Canada IT Italy Central African Republic JP Japan Congo KP Democratic People's Republic of Korea Côte d'Ivoire KR Republic of Korea Canteron LI Liechtenstein Czecheslovakia LK Sri Lanka Germany LU Laxembourg Demoark M1C Monaco	Australia FR France MN Burbados GA Gabon MR Belgium GB United Kingdom MW Burkina Faso GN Guinea NI. Bulgaria GR Greece NO Benin HU Hongary PL Brazil IE Ireland RO Canada IT Italy RU Cantral African Republic JP Japan SD Congo KP Democratic People's Republic SE Switzerland of Korea SN Côte d'Ivoire KR Republic of Korea SU Caneroon LI Liechtenstein TD Czechoslovakia LK Sri Lanka TG Germany LU Luxembourg Denmark MC Monaco

DESCRIPTION

Methods for Selecting Antibody Reagents; Anti-Idiotype Antibodies; and AIDS Vaccine Formulations

Field of the Invention

The present invention relates generally to the field of immunology and, more particularly, to novel methods of selecting for therapeutically useful antibodies. Further, the present invention relates to novel anti-idiotype anti-bodies and vaccine formulations thereof that are useful in the immunotherapy and prevention of HIV infection.

Background of the Invention

It is widely accepted that the Human Immunodeficiency 10 Virus (HIV) is the principal etiological agent of the acquired immunodeficiency syndrome (AIDS). From a clinical standpoint, depletion of CD4+ lymphocytes can account for most of the severe immunologic abnormalities in AIDS patients (Fauci, et al. (1984) Am. Int. Med. 100:92). The 15 tropism of HIV is dependent primarily upon a specific interaction between the gp120 envelope glycoprotein of HIV and the CD4 receptor molecule on target cells, e.g., T lymphocytes and macrophages, although it has been suggested that other molecules on some cells are involved in 20 similar interactions with HIV. Beyond the primary interaction between CD4 and qp120, the following steps are postulated to be essential for HIV entry into target anchorage of the virus to the cell membrane; fusion with the cell membrane; and penetration of the 25 viral elements into the cell. Any reagent that interferes with these steps is considered to block the virus infection and is of potential therapeutic value. Such reagents under current investigation include antibodies, which are generally functionally described "neutralizing 30 antibodies."

Amino acid sequence analysis of various HIV isolates from a single patient as well as different patients indicates that gp120 contains hypervariable regions inter-This finding spersed among highly conserved regions. indicates that HIV variation in vivo can occur rapidly The degree of diversity during the chronic infection. among different strains of HIV within the gp120 molecule may be due in part to natural selection during interactions with the host immune system. Antibodies that bind 10 to hypervariable regions of the HIV envelope generally only react with a single isolate or strain of the virus: such antibodies are termed "type-specific." At present, about 20 different strains of HIV have been identified. For example, antibodies directed against one of these 15 hypervariable regions of gp120 known as the "V3 loop" will neutralize HIV. However, this effect has generally been type-specific; i.e., it is well documented that among anti-gp120 antibodies, those specific to the V3 loop may exert a neutralizing effect on only a single isolate or 20 strain of HIV. Moreover, a single amino acid mutation in the V3 region allows the virus to escape neutralization by Although the precise role of the V3 such antibodies. domain in virus infection is still in question, several reports have demonstrated that antibodies to this region 25 neutralize virus infectivity by interfering with a postbinding event (Skinner, et al. (1988) J. Virol. 62:4194). Whatever the V3 region's role may be in virus infection, its variability from strain to strain and high mutation frequency make it an uncertain target for immunotherapy.

Despite the considerable sequence variation and hypervariability within gp120, all known HIV isolates share a common functional characteristic -- the ability to bind to and infect CD4+ cells. This observation leads to a basic premise that the CD4 attachment site on gp120 35 includes one or more conserved domains. antibodies binding to this CD4 attachment site could be expected to prevent viral infections and to interact fully

with various HIV isolates via their ability to bind to conserved regions of qp120. Antibodies that bind to conserved regions of HIV and recognize multiple strains of the virus are termed "group-specific." Thus, antibodies 5 that will neutralize multiple strains of HIV are said to be capable of group-specific neutralization. neutralizing antibodies would play an important role in preventing viral infection (Weiss, et al. (1985) Nature as well as killing HIV-infected cells via 10 antibody-dependent cell-mediated cytotoxicity (Rook, et al. (1987) J. Immunol. 138:1064) and complement-mediated It has been well documented that antibodies cytolysis. specific to the conserved region exist in the sera of HIVinfected humans (Schmittman, et al. (1988) J. Immunol. 15 141:4181; Moore, J.P. (1990) AIDS 4:297). Therefore, an HIV vaccine which induces antibodies capable of groupspecific neutralization may be effective in producing a clinical improvement in AIDS patients and preventing infection in healthy people.

Although the clinical correlation between neutral-20 izing antibody titer and disease progression has not been well defined, recent reports suggest that the presence of maternal antibodies to gp120 correlates positively with the uninfected status of children born to seropositive 25 mothers (Rossi, et al. (1989) Proc. Natl. Acad. Sci. 86:8055; Devash, et al. (1990) Proc. Natl. Acad. Sci. However, as most HIV-seropositive individuals have a high titer of anti-HIV antibodies as well as demonstrable levels of neutralizing antibodies, it remains 30 unresolved as to why so many of these people have a rapidly progressive infection. In attempting to answer this question, Applicant has noted in the literature that high titers of type-specific neutralizing antibody in HIVinfected individuals are not sufficient to neutralize 35 continuously emerging new isolates brought about by mutations of gp120 (Nara, et al. (1990) In Vaccines 90:137, Cold Spring Harbor Laboratory, Cold Spring Harbor, NY).

35

On the other hand, titers of group-specific neutralizing antibodies in HIV-infected individuals, which are specific against conserved regions, are low. It is, therefore, Applicant's premise that HIV vaccines must be evaluated 5 with respect to their ability to induce high titers of group-specific neutralizing antibodies.

Although the need for a HIV vaccine is clear and the objective rationally conceived, the vaccine strategies adopted have been diverse and, to date, unsuccessful. 10 Most HIV vaccine development strategies have been based on preparations including whole or portions of the major Native, recombinant, synthetic, envelope protein. vaccinia-encoded envelope proteins and whole inactivated virus have been prepared and tested.

For example, chimpanzees were immunized with the HIV encoded in a vaccinia virus. The envelope protein animals produced both humoral (antibody) and cellular (T-cell) immune responses specific for the envelope protein; however, these primates were not protected from 20 infectious challenge with HIV (Hu, et al. (1989) Nature 328:721).

Kennedy and colleagues, using synthetic peptides which represent conserved immunodominant epitopes of gp41, immunized chimpanzees with a peptide carrier conjugate 25 made from the synthetic peptides (Kennedy, et al. (1989) In Vaccines 250, Cold Spring Harbor, NY). Subsequent challenge with HIV, however, showed that the animals were not protected against infection. The use of killed whole virus vaccine (Salk, J. (1987) Nature 327:473) may be able 30 to induce protection, but to the extent these vaccines contain envelope, this protection is expected to be only type-specific and does not circumvent the problem with the envelope polymorphism (Hahn, et al. Science (1986) 23:1548).

According to a recent report (Berman, et al. (1990) Nature 345:622), immunization of primates with recombinant gp120, induced neutralizing antibodies and elicited

WO 92/15885 PCT/US92/01749

5

pr tective immunity against homologous strains of HIV. Still, there was no evidence that such immunization could provide protective immunity against different HIV strains or even mutated homologous isolates.

5

30

Antibody based vaccine approaches have also been For example, the use of anti-idiotype investigated. antibodies ("anti-ids") as vaccines against various infectious agents ("antigens") (so-called idiotype vaccines) is well documented. To date, most idiotype vaccine develop-10 ment strategies have been based on the following prepara-First, to generate monoclonal tions and principles. anti-ids (Ab2s), animal Ab1s (in most cases antibodies from non-human sources have been used) have been utilized as immunogens or templates. Second, monoclonal Abls have 15 been used as templates in the majority of cases. the in vitro selection of Ab2s as vaccine candidates was driven by the concept of antigen mimicry by anti-ids (Ab2s); i.e.. the principle of blocking the interaction of Abl with Ab2 by antigen has been used as the controlling 20 criteria for selecting Ab2 vaccine candidates. Such antiidiotype antibodies have become know in the art as "Ab2 β " or "internal image" Ab2s. Two other types of Ab2s generated against Abls have been designated: Ab2 λ , which block antigen binding to Ab1, but are not internal images 25 of the antigen; and Ab2 α , which bind to the Ab1, but do not block binding of the antigen to Abl. See FIGURE 1. Fourth, to further define Ab2s as vaccine candidates in vivo, non-primate animals have generally been used as animal models to evaluate the Ab3 response.

Despite the tremendous financial and human resources that have been dedicated to the discovery of an effective AIDS vaccine, the results have been inconclusive at best. The reasons for these apparent failures are unclear, but some of these studies indicate that broadly reactive anti-35 HIV neutralizing antibodies were not generated in sufficient titer in the experimental animals, allowing HIV to "escape" the immune response.

Summary of the Invention

The present invention is based on Applicant's discovery of a novel method for selecting antibody reagents for vaccine formulations that are capable of inducing group-neutralizing anti-HIV antibody responses in humans.

In its broadest embodiment, the present invention is directed to a method of selecting monoclonal anti-idiotype antibody reagents useful in vaccine formulations for the The first step treatment or prevention of HIV infection. in this method requires preparing human polyclonal antigp120 antibodies (Abls) and subsequently generating a first group of monoclonal anti-idiotype antibodies (G1-Ab2s) capable of reacting with or immunologically binding with the human polyclonal anti-gp120 antibodies. from this first group of monoclonal anti-idiotype antibodies (G1-Ab2s), a second group or subset of monoclonal anti-idiotype antibodies (G2-Ab2s) are selected that are characterized by their ability to immunologically bind gp120 antibodies that neutralize in vitro multiple 20 strains of HIV. Finally, the second group of monoclonal anti-idiotype antibodies (G2-Ab2s) is screened for a subset or a third group of monoclonal anti-idiotype antibodies (G3-Ab2s) which are characterized by their ability to generate an anti-anti-idiotype antibody (Ab3) 25 response in a primate host, wherein the anti-anti-idiotype antibody (Ab3) immunologically binds with gp120 and neutralizes in vitro multiple strains of HIV.

In a preferred embodiment of the invention, the Abls are prepared from pooled human sera from HIV-positive asymptomatic individuals and purified by passing over a gp120 column to increase the concentration of anti-gp120 antibodies. Also preferred is the generation of murine monoclonal Ab2s.

The present invention is further directed broadly to monoclonal anti-idiotype antibodies that are capable of eliciting a group neutralizing anti-HIV antibody response in humans. The present invention is also directed to

monoclonal anti-idiotyp antibodies (G3-Ab2s) selected acc rding to the meth ds described her in.

The present invention is further directed to vaccine formulations including the monoclonal anti-idiotype antibodies (G3-Ab2s) disclosed herein.

Thus, it is an object of the present invention to develop a method for selecting antibody reagents useful in vaccine formulations which are capable of eliciting in vivo high titers of group-specific neutralizing antibodies against HIV in humans.

It is a further object of the present invention to develop antibody reagents per se that are useful in vaccine formulations and the cell lines capable of producing such antibodies.

15 Brief Description of the Drawings

FIGURE 1 graphically illustrates the concepts of "idiotype matching" (1A) and "idiotype mimicry" (1B).

FIGURE 2 shows that the flowthrough fraction of Total Anti-gp120 Abs bound to $gp120_{IIIB}$ and not to $gp120_{SF2}$.

FIGURE 3 shows that the Total Anti-gp120 Abs bound to $gp120_{IIIR}$ and not to $gp120_{SF2}$.

FIGURE 4 graphically depicts the purification of Total Anti-gp120 Abs.

FIGURE 5 illustrates the effective inhibition of $^{125}\text{I-labeled gp120}_{\text{SF2}}/\text{ CD4}$ by CD4-site antibodies.

FIGURE 6 depicts HIV virus neutralizing activity of CD4-site antibodies.

FIGURE 7 depicts HIV virus neutralizing activity of non-V3 antibodies.

FIGURE 8 lists the first group of anti-idiotype antibodies and the characteristics of corresponding Ab1s.

FIGURE 9 summarizes the characterization of corresponding Ab3s against candidate Ab2 vaccines in monk ys.

FIGURE 10 summarizes the neutralizing activity of Ab3 isolated from 3C9 immunized monkeys.

FIGURE 11 depicts neutralizing activity of human 3C9+Ab against HIV laboratory strains (<u>Left</u>) and patient primary HIV isolates (<u>Right</u>).

FIGURE 12 depicts the binding of gp120_{SF2} (<u>Left</u>) and 5 gp120_{IIIB} (<u>Right</u>) to purified Ab3 from monkeys PR0703 and PR0783.

Detailed Description of the Preferred Embodiments

The present invention is based on Applicant's novel anti-idiotype (anti-id or Ab2) approach to developing an HIV vaccine that possesses the critical biological property of eliciting group-specific neutralizing antibodies. Although the Applicant does not wish to be bound by any particular theory, the present invention is driven by the concept of "idiotype matching" rather than the prior art concept of "antigen mimicry" by anti-ids.

"Idiotype matching" utilizes anti-ids as surrogates for antigen, e.g., in an attempt to stimulate specific elements of the B cell repertoire in humans. In contrast to the prior art, which focus exclusively on idiotype mimicry and Ab2B, the anti-id vaccine candidates of the present invention may be Ab2 A, B, or λ as long as they elicit group neutralizing Ab3 responses. Compare, FIGURE 1A and FIGURE 1B.

In short, most prior art anti-idiotype vaccine 25 development strategies have been based upon antigen mimicry and the resulting principles:

- (1) To generate monoclonal anti-ids, non-human Abls have been utilized as templates.
- (2) Monoclonal Abls have been used as templates in 30 the majority of cases.
 - (3) The <u>in vitro</u> assays used to select Ab2s as vaccine candidates have generally been based on the concept of antigen mimicry by anti-ids (Ab2B).
- (4) To further d fine Ab2s as vaccine candidates in vivo, non-primates have been the animal model of choice.

Notwithstanding th scientific momentum created by these teachings in the art, Applicant develop d an entirely new approach for the selection of anti-id vaccines. Some of the preferred features and advantages of 5 Applicant's invention include the following:

- (1) Human Abls, preferably purified from seropositive healthy individuals, are used as templates to generate Ab2s. This choice of starting material increases the likelihood that the Ab2-based vaccine will generate 10 Ab3s in humans that are immunologically equivalent to the human Abls with respect to antigen binding potential, This feature of the invention i.e., anti-HIV binding. takes advantage of the fact that Applicant knows that B cell repertoires pre-exists in humans that are capable of 15 producing Ab1-like Ab3s (so called Ab1') In contrast, using a murine Ab1 assumes that the B cell repertoire of the mouse is genetically and immunologically equivalent to the human B cell repertoire--possibly a fatal flaw for vaccine development.
- (2) Polyclonal Abis rather than monoclonal Abis are used as templates to generate Ab2s. This feature of the invention eliminates the possibility that an individual monoclonal Ab1 (which represents only one epitope) selected for use as a templates does not represent a 25 highly conserved HIV epitope.
- Applicant has discovered that the definition of Ab2 B and the evaluation of Ab3 responses in non-primate animals are not the critical factors in selection of an Ab2 as a vaccine candidate. In contrast, Applicant has 30 discerned the following criteria for the selection of anti-idiotype vaccines: (a) Abls (which react with the Ab2 vaccine candidate) from the human sera exhibit desirable biological effects in vitro (e.g., virus neutralization). Ergo, Ab3s generated in vivo in 35 primates, including humans, in response to an Ab2 vaccine are expected to have Ab1-like biological effects. Ab2s induce group-specific neutralizing Ab3s in primates.

The principal advantages in using Applicant's anti-id approach to select an HIV vaccine are threefold. First, selected anti-ids elicit only beneficial anti-gp120 anti-bodies (e.g., group-specific neutralizing antibodies).

5 Second, anti-ids elicit high titers of certain anti-gp120 antibodies. Third, the biological properties of anti-gp120 antibodies elicited by these Ab2 are predictable.

in its broadest embodiment, the present invention is directed to a method of selecting monoclonal 10 anti-idiotype antibody reagents that are useful in vaccine formulations for the treatment and prevention of HIV The first step in this method involves the infection. preparation of human polyclonal anti-gp120 antibodies (Abls) that are to be used as templates to generate Ab2s. 15 It is preferred that the human polyclonal anti-gp120 antibodies are obtained from human sera of HIV-positive, asymptomatic individuals. It is further preferred that such sera is purified to increase the representation of anti-gp120 antibodies by passing the pooled human sera 20 over a gp120 column. It is also contemplated by the present invention that the total anti-gp120 antibodies will be passed over not one but two gp120 columns, wherein each gp120 column represents a different strain of gp120. It is most preferred that the human polyclonal anti-gp120 25 antibodies represent the total anti-gp120 antibodies. Alternatively, the total anti-gp120 antibody population can be dissected into fractions for use in the present attachment site-specific anti-gp120 CD4 invention: antibodies, V3 region-specific anti-gp120 antibodies or 30 non-V3 region-specific anti-gp120 antibodies.

The second step of the method of the present invention involves the generation of a first group of monoclonal anti-idiotype antibodies which are capable of reacting with the human polyclonal anti-gp120 antibodies pr pared in the first step of the method. Preferably, these anti-idiotype antibodies are murine monoclonal antibodies generated by conventional hybridoma technology.

As used herein, the term "monoclonal" is intended to enc mpass antibodies derived from a singl clone as well as mono-specific antibodies derived by molecular biology techniques. Alternative sources of anti-idiotype antibodies contemplated by the present invention include human, primate, chimeric, and humanized antibodies, as well as antibodies synthesized by molecular biology techniques. As used herein, the term "antibody" is intended to include whole antibodies, antibody fragments, 10 Fab fragments, single chain antibodies, and immunologically functional antibody equivalents.

Once the first group of monoclonal anti-idiotype antibodies have been generated, the third step in the method of the present invention involves screening this 15 first group of antibodies for a second group (subset of first group) of monoclonal anti-idiotype antibodies that are characterized by their ability to react with or immunologically bind anti-gp120 antibodies that neutralize In a preferred in vitro multiple strains of HIV. 20 embodiment, the anti-gp120 HIV neutralizing antibodies are human. In a most preferred embodiment, these human antigp120 HIV neutralizing antibodies are human antibodies selected from the pool of human polyclonal antibodies prepared according to the first step of the method. 25 used herein, neutralization of more than two strains of HIV is considered to be "multiple" strain neutralizing. These anti-gp120 antibodies may be selected by passing the human polyclonal antibody preparation made in the first step of the method over an Ab2 (group 1) column, and 30 selecting those anti-gp120 antibodies that are capable of in vitro neutralization of multiple strains of HIV.

The present invention is also directed to antibodies, in particular, to monoclonal anti-idiotype antibodies that are capable of eliciting a group- neutralizing anti-HIV antibody response in humans. Preferred are monoclonal anti-idiotype antibodies selected according to the selection method described above.

SUBSTITUTE SHEET

35

Hybridoma cell lines are also an aspect of the In particular, hybridoma cell lines present invention. capable of producing a monoclonal anti-idiotype antibody that is capable of eliciting a group neutralizing anti-HIV 5 antibody response in humans. The most preferred is the hybridoma cell line deposited with the ATCC, Rockville, MD on March 8, 1991, designated accession no. HB10701. addition to the antibodies produced by this hybridoma, the present invention includes monoclonal antibodies capable 10 of eliciting a substantially equivalent immune response in humans as a monoclonal antibody produced by the hybridoma cell lines described above. Preferably, such a monoclonal antibody would react with and compete for binding to the same anti-gp120 antibody as a monoclonal antibody produced 15 by the hybridoma cell lines described above.

An additional embodiment of the invention is the development of vaccine formulations including, in mixture with a pharmaceutically acceptable carrier, a monoclonal anti-idiotype antibody that is capable of eliciting a 20 group neutralizing anti-HIV antibody response in humans. Vaccine formulation, including selection of appropriate adjuvants and carriers, is well known to those of skill in the art. Once formulated, a vaccine's dosimetry would be based upon patient and clinical factors that are well 25 appreciated by those of skill in the art of vaccine deliv-The vaccines of the present invention may be either prophylactic or therapeutic in nature. Although a single anti-idiotype antibody reagent vaccine is desirable, it is contemplated by the present invention that multiple anti-30 idiotype antibodies may be formulated in a "cocktail" vaccine approach. Again, clinical factors will dictate the vaccine formulation and dosage. However, a preferred vaccine formulation would include about 2.5 mg of an antibody of the present invention formulated with about 100 ug of termutide and SAF (Syntex adjuvant formulation -- 2.0 mg polysorbate 80, 12.5 mg Poloxymer 401, and 50 mg

20

squalane). An alternative formulation may exclude the termutide component.

EXAMPLE 1A Preparation of Starting Materials -- Human Polyclonal Total Anti-gp120 Antibodies.

Total anti-gp120 antibodies (anti-gp120_{sp2}) ("Total Anti-gp120 Abs") were purified from human sera by affinity $gp120_{SF2}$ -Sepharose chromatography using a Recombinant $gp120_{SF2}$ was derived from engineered CHO cells. A pool of sera from HIV-infected asymptomatic individuals 10 (purchased from North American Biologicals, Inc., Miami, Florida 33169) were inactivated by detergent (1% NP40) and heat (1 hr. at 56°C) treatment. The sera was dialyzed against PBS and ammonium sulfate (40% as final concentration) was added to the sera to isolate the immunoglobulin 15 fraction. The precipitated immunoglobulin fraction was reconstituted with PBS, filtered and applied to the gp120_{SF2}-Sepharose column. The bound Total Anti-gp120 Abs on the column were eluted with pH 3.0 glycine buffer and dialyzed against PBS. The recovery of total anti-gp120 antibodies in sera was approximately 300 μ g/ml.

Antibodies binding to gp120_{SF2} are expected to consist of two kinds: (a) SF2 type-specific antibodies, and (b) antibodies binding to conserved regions of gp120 such as the CD4-attachment site. Since the individuals contributing to the Ig pool might have been infected with a variety of HIV strains, the flowthrough fraction may bind to gp120 derived from strains other than SF2. Indeed, as shown in FIGURE 2, Applicant found that the flowthrough fraction bound to gp120_{IIIB}, not to gp120_{SF2}, indicating the presence of III_B type-specific anti-gp120 antibodies in the flowthrough fraction. In addition, Total Anti-gp120 Abs bound to both gp120_{SF2} and gp120_{IIIB} (FIGURE 3), indicating the presence of a mixture of type-specific and cross-reactive anti-gp120 antibodies.

Applicant then investigated the neutralizing activi-35 ties of Total Anti-gp120 Abs against various HIV isolates.

FIGURE 2 shows that Total Anti-gp120 Abs exhibited neutralizing activities against four different HIV isolates, although their neutralizing activities against the different HIV isolates varied as much as 10-50 fold, indicating that Total Anti-gp120 Abs contains significantly different populations and amounts of type- and group-specific neutralizing antibodies. From these results, Applicant reasoned that the difference in the neutralizing activities of Total Anti-gp120 Abs against the various viral isolates could be explained by (a) strain-to-strain differences in the amount of anti-gp120 antibodies required for neutralization of each virion, and (b) the presence of type- or group-specific neutralizing antibodies.

To demonstrate the presence of type-specific neutralizing antibodies in Total Anti-gp120 Abs, Applicant tested the binding of Total Anti-gp120 Abs to type-specific neutralizing epitopes derived from the V3 region of different gp120s.

20 <u>Table 1</u>
<u>Binding Activities of Different Anti-gp120 Antibody</u>
<u>Preparations To Peptides Derived From the V3 Region of</u>
Various <u>qp120s</u>

				Binding A	ctivities*	
25	Anti-gp120 Antibodies			Origin of	Peptides	
		· MN	IIIB	RF	sc	WMJ-2
	Total Anti- gp120 Abs	>3	1.05	0.06	0.8	1.48
30	CD4-site Abs	0.15	0.01	0.01	0.03	0.01

*Binding activities are expressed by O.D. obtained at $5\mu g/ml$ of each preparation.

Table 1 shows that Total Anti-gp120 Abs bind extensively to the MN isolate derived epitope, extensively to the IIIB, SC and WMJ-2 derived epitopes and Thus, the highest weakly to the RF-derived epitope. 5 binding activities of Total Anti-gp120 Abs against the MNderived epitope is consistent with the fact that Total Anti-gp120 Abs exhibited the highest neutralizing activities against the MN isolate. However, it is interesting to note that the binding activity of Total Anti-gp120 Abs 10 to the RF-derived epitope was only slightly above background, despite the fact that Total Anti-gp120 Abs effectively neutralized the RF isolate. This could be explained by the presence of group-specific neutralizing antibodies specific against non-V3 epitopes.

15 EXAMPLE 1B Preparation of CD4 Site Attachment Site-Specific Polyclonal Anti-gp120 Antibodies from Human Sera

<u>Isolation and Characterization of CD4 Attachment Site-</u> <u>Specific Anti-qp120 Antibodies</u>

Since Applicant observed group-specific neutralizing 20 activities with Total Anti-gp120 Abs, it was of interest to determine the epitope specificity of anti-gp120 antibodies which are responsible for these group-specific neutralizing activities. Since antibodies to the conserved CD4 attachment site of gp120 are considered to be 25 the major group-specific neutralizing epitope, Applicant devised the following strategy. The approach to isolate CD4 attachment site-specific anti-gp120 antibodies (CD4site Abs) from Total Anti-gp120 Abs consisted of saturating the CD4-attachment site on gp120 immunoabsorbent with 30 sCD4 and followed by applying Total Anti-gp120 Abs to the immunoabsorbent in the presence of excess sCD4. flowthrough fraction of the antibodies were collected and applied to a protein G column to isolate the IgG fraction from sCD4. This is graphically depicted in FIGURE 4. To 35 ensure complete r moval of the sCD4 from the IgG fraction,

Applicant analyzed this preparation by SDS gel electrophoresis and ELISA, using a monoclonal anti-CD4 antibody, and could not detect any sCD4 by either method.

To isolate the CD4 attachment site-specific antigp120 antibodies, 10 mg of sCD4 in PBS was repeatedly
loaded on the gp120_{SF2}-Sepharose column to saturate the CD4
attachment site onto the gp120. Then, 4 mg of Total Antigp120 Abs with 10 mg of sCD4 in PBS was loaded onto the
column. The flowthrough fraction was further separated on
a protein G column (Pharmacia Fine Chemicals, Piscataway,
NJ) to isolate the IgG fraction. SDS gel electrophoresis
and ELISA, using a monoclonal anti-CD4 antibody, were
employed to verify that the isolated IgG was free of sCD4.

To examine whether the IgG fraction contains CD4-site Abs, Applicant performed competitive inhibition assays on CD4+ cells using \$^{125}I\$-labeled gp120_SF2\$ and CD4-site Abs. FIGURE 5 shows that CD4-site Abs effectively inhibited binding of \$^{125}I\$-labeled gp120_SF2\$ to CD4 cells. The approximate IC50 of CD4-site Abs in this assay was ten times less than that of sCD4. Applicant observed almost the same results with \$^{125}I\$-labeled gp120III_B. This result indicated that this IgG fraction contained significant amounts of CD4-site Abs. To determine what proportion of this IgG fraction are CD4-site Abs, various anti-gp120 antibodies were coated on microtiter plates and incubated with \$^{125}I\$-labeled gp120_SF2\$ in the presence or absence of sCD4.

17

Table 2 Proportion of CD4 Attachment Site-Specific Antibodies in Different Anti-qp120 Antibody Preparations

	Anti-gp120	Proportion*
5	Antibodies	as %
	Total Anti-gp120 Abs	16
	CD4-site Abs	71
	Non-V3 Abs	32
	V3-Region Abs	0

10 *The percent inhibition by sCD4 in each anti-gp120 antibody binding to gp120 was calculated. The percent inhibition was considered to be the proportion of .CD4 attachment site-specific anti-gp120 antibodies.

Table 2 shows that sCD4 inhibited the binding of CD4-15 site Abs to $gp120_{SF2}$ by 71%, whereas only 16% of Total Anti-gp120 Abs binding to gp120_{SF2} was inhibited by sCD4. This result indicated that our purification significantly enriched CD4-attachment site-specific anti-gp120 antibodies; however, these antibodies still contained 20 substantial amounts of non-CD4 attachment site-specific anti-gp120 antibodies. In an attempt to determine the specificities of the contaminating antibodies in CD4 site Abs, Applicant tested the binding of CD4 site Abs to synthetic peptides derived from the V3 domain of different 25 HIV isolate gp120s. Table 1 shows that CD4 site Abs bound only weakly to the MN isolate-derived epitope and not to the other isolate-derived epitopes. This finding demonstrates that the purification for CD4-site Abs removes most of the V3 region-specific anti-gp120 antibodies.

Since Applicant found that CD4-site Abs contain mainly CD4-attachment site-specific anti-gp120 antibodies with minimal amounts of MN isolate specific neutralizing antibodies, Applicant investigated whether CD4-site Abs could neutralize different HIV isolates. As shown in 35 FIGURE 6, CD4-site Abs exhibited almost the same extent of neutralizing activities against three different HIV isolates. Compared with CD4 site Abs, Total Anti-gp120

Abs exhibit much higher neutralizing activities against the MN isolate, a little higher against the RF isolate and almost the same against the ${\rm III}_{\rm B}$ isolate. From the fact that CD4 site Abs effectively neutralized three different 5 HIV isolates at or near equivalent antibody concentrations despite low or no binding activities against identical V3 specific peptides, Applicant concluded that CD4-site Abs exhibit mainly non-V3 group-specific neutralizing activities, whereas, Total Anti-gp120 Abs exhibit type- and 10 group-specific neutralizing activities. Applicant utilized these CD4-site Abs to generate murine monoclonal anti-ids.

EXAMPLE 1C. Preparation of Anti-qp120 Antibodies Depleted of Type-Specific Neutralizing Antibodies

Another approach for identifying group-specific 15 neutralizing antibodies was to deplete type-specific neutralizing antibodies from TAGA, using a column of V3 peptide Affi-Gel 401 and test this preparation for groupspecific neutralizing activities.

To separate the anti-gp120 antibodies depleted of V3 region-specific antibodies (non-V3 Abs) from Total Antigpl20 Abs, a mixture of two peptides (8 mg of each) derived from the V3 region of gp120_{SF2} (see peptides for sequence) were coupled to 10 ml of Affi-Gel 401 (Bio-Rad, 25 Richmond, CA). See method below. The peptides were synthesized on an automated peptide synthesizer (model 430A, Applied Biosystems, Foster City, CA) by solid phase The peptides were deprotected and cleaved from method. the support with HF/anisole or the "low-high" HF proce-The sequences of two SF2-derived peptides for an 30 dure. **V3SF2:** follows; were as column affinity IGDIRKAHC and TRPNNNTRKSIYIGPGRAFHTTGRI The sequences of CTRPNNNTRKSIYIGPGRAFHTTGRIIGDIRKAHC. region peptides for binding assays w re as various V3 YNKRKRIHIGPGRAFYTTKNIIG, <u>MN</u>: 35 follows; NNTRKSITKGPGRVIYATGQIIG, NNTRKSIRIORGPGRAFVTIGKIG, RF:

SC: NNTTRSIHIGPGRAFYATGDIIG, WMJ-2: NNVRRSLSIGPGRAFRTREGN. Ten mg of Total Anti-gp120 Abs in PBS was repeatedly loaded on the peptide Affi-Gel column. The flowthrough, consisting of anti-gp120 antibodies depleted of V3 region-specific antibodies (non-V3 Abs), was collected. SF2 V3 region-specific anti-gp120 antibodies (V3-region Abs) were eluted with pH3 glycine buffer.

Peptide Coupling

10 1. Derivatization of the Affinity Support

An aliquot of Affi-Gel 401 affinity support was reduced with 100 mM dithiothreitol for 1 hour at room temperature. The resin was then transferred to a sintered glass funnel and washed extensively with argon-saturated phosphate buffered saline.

The reduced resin was then transferred to a 15 ml conical tube and diluted with 10 ml of argon-saturated phosphate buffered saline. The bis-maleimidomethyl ether in ethanol was added to the resin in a 100 fold molar excess of maleimide over sulfhydryls. The reaction was allowed to proceed overnight at 1-8°C with rotation in a sealed tube. The resin was subsequently transferred to a sintered glass funnel and washed extensively with argon-saturated PBS in order to remove excess reagent.

25 Unreacted sulfhydryls were blocked by the addition of n-ethylmaleimide in ethanol in a sealed tube with rotation at 2-8°C overnight.

2. Reduction of Peptide and Coupling to the Derivatized Support

Reduction of the C terminal cysteine was effected by the addition of 100 mM DDT with incubation for 2 hours at room temperature. The reduced peptide was separated from r ductant by adsorption to a Sep-Pak C-18 apparatus, washed with argon-saturated PBS, and eluted directly into the support in solution with 95% CH₃CN. Coupling was

allowed to take place overnight at 2-8°C with rotati n in a sealed tube. Following the coupling, the resin was washed with PBS in order to remove excess peptide.

Both bound and flowthrough fractions were tested in 5 virus neutralization assays. FIGURE 7 shows that non-V3 Abs exhibited almost the same extent of neutralizing activities against three different HIV isolates. However, V3 region Abs exhibited strong neutralization activities against the MN isolate, minimal activities against III_B and 10 no activity against RF. Based on these observations, Applicant concluded that the neutralizing activities of non-V3 Abs are group-specific while the neutralizing activities of V3 region Abs are type-specific. result indicates that the majority of neutralizing 15 activity of Total Anti-gp120 Abs to III_B and RF isolates (FIGURES 7a, 7b) are derived from group-specific neutralizing activities of the non- V3 Abs fraction, whereas, most of the neutralizing activity of Total Anti-gp120 Abs to the MN isolates (FIGURE 7c) are derived from type-20 specific neutralizing antibodies in the V3 region Abs Furthermore, by comparison of neutralizing fraction. activities of V3 region Abs and non-V3 Abs on MN isolate (FIGURE 7c), Applicant concluded that V3 region Abs are more effective than non-V3 Abs neutralizing antibodies in 25 neutralizing a specific isolate.

Conclusions on Preparation of Anti-gp120 Antibodies in Examples 1A, 1B, and 1A

Applicant has dissected the neutralizing activities of various polyclonal anti-gp120 antibodies present in a pool of serum from 4 HIV-infected individuals to understand the role of these antibodies in the host defense against HIV infection. The major findings of this study are as follows: (i) Total Anti-gp120 Abs exert type- and group-specific neutralizing activities; (ii) the group-specific neutralizing activity is found primarily in CD4-site Abs; (iii) if other group-specific neutralizing

10

20

antibodies rather than CD4- site Abs exist, they are not specific against a simple linear epitope f the V3 domain; (iv) V3 region Abs are more effective than non- V3 Abs in neutralizing a specific HIV isolate.

Although neutralizing efficacy of CD4-site Abs is relatively low, it is reasonable to assume that an HIV vaccine which induces high titers of CD4 site Abs should be effective in producing a clinical improvement in AIDS patients and preventing infection against various HIV isolates in healthy people.

EXAMPLE 2A Generation of First Group of Anti-id Mabs and Characterization for their Ability to React with Human Polyclonal Anti-qp120 Antibodies.

Generation of anti-id Mabs was accomplished by fusing 15 sp2 murine myeloma cells with splenocytes from mice which were immunized 2-3 times with polyclonal anti-gp120 antibodies. In screening assays, anti-id was selected for (i) binding to polyclonal human anti-gp120 antibodies, (ii) not binding to any specific isotype standard antibody, (iii) no binding to Ig preparation of the same pool of sera depleted of Total Anti-gp120 Abs by absorption on a gp120_{SF2}-Sepharose column.

<u>Immunizations</u>

BALB/c mice were immunized with TAGA and the various 25 other fractions of polyclonal anti-gp120 antibodies purified from a pool of four sera from HIV-infected individuals as described in EXAMPLES 1A-C. All immunizations were done by two S.C. injections of 50 micrograms of antibody in SAF-MDP adjuvant at two week intervals, followed 30 by a single i.p. injection of 50 micrograms of antibody in incomplete Freunds adjuvant ("IFA") 3 days prior to Spleen cells were fused with SP2 myeloma cells and hybrids were selected in HAT medium. Culture supernatant were tested by ELISA. Positive clones were

subclon d twice. Anti-id Mabs were purified from the ascites on Protein A-Sepharose columns.

ELISA Assay for Hybridoma Screening

For hybridoma screening assays, microtiter wells were coated with 300 ng of the following human Igs in PBS: Different anti-gp120 antibodies used for immunizations; Ig preparations of the same pool of sera depleted of Total Anti-gp120 Abs by absorption on a gp120_{SF2}-Sepharose column; and myeloma proteins representing all different isotypes of human heavy and light chains (The Binding Site, San Diego, CA). Culture supernatants were tested for their binding activities to the antibody coated plates. Alkaline phosphatase-coupled goat anti-mouse Ig and phosphatase substrate (Sigma, St. Louis, MO) were used to detect antibodies binding to coated plates.

EXAMPLE 3A. Selection of a Second Group of Anti-IDS from the First Group of Anti-IDS Based on in vitro Characterization

Characterization of anti-id Mabs was necessary to identify appropriate antibodies for the next level of experiments. For this procedure, a unique approach to define Ab2s as vaccine candidates in vitro was employed. The underlying principals are based on "idiotype matching" rather than the antigenic mimicry of anti-id (Ab2B).

Therefore, the isolation of the corresponding Ab1s to individual Ab2s from human sera and their evaluation was a crucial step. Thus, Applicant isolated the corresponding Ab1s using an Ab2-Sepharose column and tested the Ab1s for binding activity to gp120 and neutralizing activity against multiple strains of HIV. Then, Ab2s which interact with group-specific neutralizing Ab1s were selected for primate studies. See FIGURE 8.

Preparation of Id+Abs (Abls)

Anti-id Mabs wr conjugated to Sepharose 4B by mixing 8 mg of antibodies with 5 ml of CNBr-activated Sepharose 4B (Pharmacia, Piscataway, NJ). Then, 20-60 ml of the same sera pool used for various anti-gp120 antibody preparation was loaded to anti-id Mab-Sepharose 4B columns. After washing extensively, bound antibodies (Abl Ids) were eluted with citrate buffer ph2.8 and dialyzed against PBS.

10 ELISA

An ELISA plate was coated overnight with 0.2 ml of $0.5~\mu \rm g/ml$ of $\rm gp120_{SF2}$ or $\rm gp120III_B$ (SmithKline Beecham, King of Prussia, PA) in 0.1 M bicarbonate buffer pH 9.8. The plate was washed three times with PBS and blocked for 15 1 hr. with 10% FCS in PBS. After washing, various concentrations of Ab1s were added. Two hours later, the plate was washed and biotin-coupled goat anti-human Ig and avidin-coupled peroxidase was used for detecting antibody binding to the ELISA plate.

20 <u>Neutralization assay</u>

A quantitative neutralization assay of HIV by antigp120 antibody was conducted as previously described (Nara, et al. (1988) Nature 332:469). The virus stocks of HIV isolates were grown in H9 cells and harvested for optimal viral infectivity and aliquoted and frozen at -140°C in a vapor phase cryo unit. Equivalent amounts of virus (syncytial-forming unit) were used in the assays so as to allow for comparison of the various antibody fraction tested. The assays are repeatable over a virus surviving fraction range of 1 to 0.001 within a two to four fold dilution of the antibody (p≤0.001).

EXAMPLE 4A. Selection of Third Group of Anti-IDS from the Second Group Based on AB3 Response in Primates

Demonstrations of HIV Vaccine Activity of Anti-Id Mabs

To identify anti-id Mabs with potential utility in eliciting group-specific neutralizing anti-gp120 anti-bodies, cynomolgus monkeys were used. Monkeys were immunized with anti-id Mab/Syntex Adjuvant Formulation ("SAF") (Syntex Corp., Palo Alto, CA) emulsion, each given at two week intervals. After the fourth immunization, blood samples were taken and sera were isolated and used to prepare Ab3 preparation on Ab2 affinity column. Applicant next examined whether Ab3 binds to gp120.

Ab3 binding to gp120

To test gp120 binding activity of Ab3, gp120_{SF2} and gp120_{IIIB} were iodinated by Bolton-Hunter reagent (ICN Biochemical, Cleveland, OH) and lactoperoxidase (Bio-Rad, Richmond, CA) respectively. Microtiter plates were coated for 18 hours with 1µg of Ab3 in PBS. After washing, the plates were blocked for one hour with 1% BSA in PBS, 10,000 cpm of gp120_{SF2} or gp120III_B were added to wells and incubated for 3 hours. After washing three times, each well was counted.

HIV Neutralization

For neutralizing assays, Ab3 was further purified on gp120-Sepharose affinity column to enrich gp120 specific Ab3s. Applicant then tested whether gp120 specific Ab3s neutralize multiple strains of HIV.

Two quantitative neutralization assays of HIV were performed, as previously described (Nara, et al. (1988) Nature 332:469). Briefly, target cells of the CEM-SS cell line were infected in monolayers by different HIV isolates grown in H9 cells, in the presence of Ab3s. After 7-12 days in culture, syncitia forming units were counted, and supernatants were assayed in ELISA for p24 concentration.

As shown in FIGURES 9 and 10, Applicant found that Ab3s isolated from the sera of 3C9 immunized monk ys bind to gp120 and the purified gp120-specific Ab3s neutralize three distinctive HIV isolates.

Characterization of 3C9 EXAMPLE 4B. Immunotherapy Reagent <u>Immunizations</u>

Cynomolgus monkeys received intramuscular Four injections with 2.5 mg of 3C9 in 1 ml of SAF (Syntex containing 0,6 of formulation) 10 adjuvant N-acetylmuramyl-L-alanyl-D-isoglutamine (MDP). The injections were administered bi-weekly. Animals were bled before the injection series and seven days after each injection.

Recombinant Proteins 15

20

gp120_{IIIB} is a recombinant glycoprotein of an HIV_{IIIB} isolate which is secreted from a Drosophilia cell line. gp120_{SF2} is a recombinant envelope glycoprotein of an HIV_{SF2} isolate which is secreted from Chinese hamster ovary The $gp120_{IIIB}$ and $gp120_{SF2}$ are readily (CHO) cells. available to those of skill in the art. Recombinant soluble CD4 (sCD4) was prepared by transfection of a truncated CD4 gene into CHO) cells. Secreted sCD4 in culture supernatants of CHO cells (representing 95% of the 25 amino-terminal extracellular portion of the molecule) was purified on an anti-CD4 antibody conjugated Sepharose column.

Affinity Purification of Antibodies

All affinity columns were prepared by conjugating: 8 mg of protein with 5 ml of CNBr-activated Sepharose 4B 30 (Pharmacia). To purify the 3C9-reactive antibody (3C9+Ab) from pooled sera of healthy HIV-infected individuals (Kang, et al. (1991) Proc. Natl. Acad. Sci. USA 88, 6171-6175), heat (56°C, 1hr)- and detergent (1% Nonidet P-40)-

treated sera were applied to a 3C9-conjugated Sepharose affinity column. After being washed extensively, bound antibodies (3C9+Abs) were eluted with citrate buffer (pH 2.8) and dialyzed against phosphate-buffered saline (PBS). To purify various Ab3s from monkey sera, all 3C9 immune sera were initially diluted 1:1 with PBS. The diluted serum was then passed over an affinity column of normal Sepharose. immunoglobulin-conjugated mouse flow-through fraction was collected and set aside. After 10 washing, the bound antibody (anti-isotype Ab) was washed, eluted with citrate buffer, and dialyzed against PBS. Subsequently, the Ab3 was passed over a gp120sr-Sepharose affinity column, and the column was washed extensively. The bound antibody (gp120-specific Ab3) was eluted with 15 citrate buffer and dialyzed against PBS.

RIA

gp120_{SF2} and gp120_{IIIB} were labeled with ¹²⁵I by using a Bolton-Hunter Reagent Kit (NEN). RIA plates (Dynatech) were coated for 18 hr with various amounts of antibody in 20 PBS. The plates were washed and then blocked for one hr with PBS containing 1% bovine serum albumin (BSA). ¹²⁵I-labeled gp120_{IIIB} or gp120_{SF2} in PRS containing 10% (vol/vol) fetal calf serum and 0.05% Tween 20 was added in the presence or absence of different concentrations of inhibitors. After 3 hr, the plates were washed with PBS and the bound radioactivity was determined in a gamma counter.

ELISA

To test the binding activity of 3C9+Abs to a V3
region-derived peptide (amino acids 301-338 of gp120_{SF2}),
microtiter plates were coated overnight with 0.25 ug of
protein in PBS per well. After the adherent layer was
blocked and washed, different concentrations of antibody
in PBS were added to the plates. Alkaline phosphatasecoupled goat anti-human immunoglobulin and phosphatase

WO 92/15885 PCT/US92/01749

27

substrat (Sigma) were us d to d tect antibodies binding to th plates.

Neutralization Assays

Two quantitative neutralization assays with HIV-1

1 aboratory stains were performed, as previously described
(Nara & Fischinger, (1988) Nature (London) 332, 469-470:
Whalley et al. (1992) Viral Immunol., 4, 201-213).
Briefly, monolayers of CEM-SS target cells were cultured with different HIV strains, in the absence or presence of
different neutralizing antibodies. After 3-5 days in culture, syncytia-forming units were counted and in the same assay, supernatants were assayed in ELISA for p24 concentration 7-12 days later. Neutralization assays with HIV primary isolates were performed as described (Daar et al. (1990) Proc. Natl. Acad. Sci. USA 87, 6574-6578).

Property of Anti-Idiotypic mAb 3C9 as Candidate for Clonotypic Stimulation

Using the methods described above, Applicant generated a series of anti-idiotypic antibodies against 20 purified human anti-gpl20 antibodies from a pool of four healthy HIV-1 infected donors. Briefly, anti-gp120 antibodies were purified on gp120_{sp2}. Sepharose and used to immunize mice for hybridoma fusion. Anti-idiotype mAbs were selected for binding to human anti-gp120 antibodies. 25 One anti-idiotype mAb, 3C9, with specificity for CD4 site-directed anti-gp120 antibodies, was characterized as a candidate for induction of anti-gp120 antibodies. Pooled human HIV-1+ sera were purified on a 3C9 immunoabsorbent and eluted antibodies (3C9+Ab) were 30 tested for binding to gp120 and V3 loop peptides. results (Table 3), show that 3C9+Ab bound to two genetically distinct gp120s derived from ${\rm HIV}_{\rm SF2}$ and ${\rm HIV}_{\rm HIB}$ and did not bind to the V3 loop peptide. Furthermore, 3C9+Ab binding to gp120 was inhibited by sCD4, indicating 35 that 3C9+Ab is specific for an epitope around the CD4

attachment site of gp120. Next, the biological function of 3C9+Ab was examined in virus neutralization assays. Human 3C9+Ab exhibited different spectra of broadly neutralizing activities against genetically distinct HIV 5 laboratory strains and four primary isolates established individuals. See Piq. 11. HIV-infected Collectively, these results indicate that anti-idiotype mAb 3C9 recognizes a clonotypic determinant shared by human polyclonal anti-gp120 antibodies that are specific 10 for a conserved epitope around the CD4 attachment site and exhibit broadly neutralizing activities. definition, 3C9 defines a marker on B cells producing neutralizing anti-gp120 antibodies.

29

Tabl 3.

Immunol gical Prop rti s of Human 3C9-Reactive
Antibodies (3C9+ Ab)

Binding activity of 3C9+ Ab in assays cpm in RIA with 125I-qp120†

A in ELISA

		qp12	<u>0</u> sf2 _	gp120 ₁₁	tw 81	th V3
10		Without	With	Without	With	loop
	Ab coating*	sCD4	sCD4	sCD4	sCD4	peptidett
	3C9+ Ab	5058	1539	8394	1963	0.07
	Normal hIg	, 145	151	282	278	
	V3-specific hAB	3		***		1.77
15	(polyclonal)					

†Plates were coated with 0.3 ug of antibody per well. After addition of PBS containing 10% fetal calf serum to block the absorbed layer, 5 x 10⁴ cpm of ¹²⁵I-labeled 20 gp120_{SF2} or ¹²⁵I-labeled gp120_{IIIB} was added to each well in the presence or absence of sCD4 at 50 ug/ml. After 3 hr, the plates were washed and the radioactivity was measured.

3C9+ Ab was purified from HIV human sera on a 3C9-conjugated Sepharose column, normal human immuno-globulin (hIg) was from The Binding Site (San Diego), and V3-specific human polyclonal antibodies (hAb) were purified from HIV human sera on a V3_{SF2}-conjugated Affi-Gel column as described (Kang et al. (1991) Proc. Natl. Acad. Sci USA 88, 6171-6175).

30 ††Plates were coated with 0.25 mg of V3_{SF2} peptide in PBS per well. After addition of PBS containing 1% BSA to block the absorbed layer, 1 ug of antibodies per ml of PBS were added to the plates. Alkaline phosphatase-coupled goat anti-human Ig and phosphatase substrate were used to detect antibodies binding to the plates. Absorbance was measured at 405 nm.

3C9 mAb Induces Anti-qp120 Ab3 Antibodies

To test whether 3C9 elicit anti-gp120 antibodies, four cynomolgus monkeys were immunized with 3C9 mAb in Syntex adjuvant formulation. Since the immune serum was 5 expected to contain 3C9-specific antibodies (Ab3) and mouse immunoglobulin-specific antibodies (anti-isotype Ab), it was fractionated on successive affinity columns and tested for binding to gp120. Ab3 from two monkeys bound to gp120_{SF2} in an antibody concentration-dependent 10 manner and to gp120IIIB in an antigen concentrationdependent manner, whereas anti-isotype Abs did not bind to either type of gp120. See Fig. 12. Similar results were obtained with Ab3s from two other monkeys (data not shown). Inhibition assays were performed to analyze the The results in Table 4 show that 15 specificity of Ab3. gp120 binding of Ab3 from monkeys PR0703 and PR0783 were inhibited by 3C9 anti-idiotype mAb (Ab2), human 3C9+Ab (Abl), sCD4. Similar results were obtained with Ab3s from other monkeys (data not shown). These results confirm the 20 anti-clonotypic nature of 3C9 as it is able to induce in primates gp120-specific Ab3 with similar specificity to the human anti-gp120 antibodies (3C9+Ab) purified on 3C9 immunoadsorbent (See Table 3).

WO 92/15885 PCT/US92/01749

<u>;</u> ·

Table 4 Inhibition of Ab3 Binding to ${\tt gp120_{SF2}}$ by Vari us Inhibitors

5	Monkey Ab3 bindir	ng		
	to qp120 _{SF2}	Inhibitor	Conc., ug/ml	%inhibition
	PRO703 Ab3	3C9	5	94
		3C9+ Ab+	5	98
		sCD4	30†	73
10		BSA	100	2
	PRO783 Ab3	3C9	5	95
		3C9+ Ab+	5	97
		sCD4	30†	79
		BSA	100	0
15				

Plates were coated with 2 ug of antibody per well. After addition of PBS containing 10% fetal calf serum to block the adsorbed layer, 1.5 x 10⁵ of ¹²⁵I-labeled gp120_{SF2} was added to each well in the presence or absence of inhibitors. After 3 hr, the plates were washed, the radioactivity was assayed, and the percent inhibition was

calculated.

*3C9⁺ Ab was purified from HIV⁺ human sera on a 3C9-conjugated Sepharose column.

25 †At this concentration, sCD4 did not inhibit V3 region-specific antibody to $gp120_{SF2}$, as previously described.

Monkey Ab3 antibodies were further fractionated on a gp120_{SF2}-Sepharose affinity column. The flowthrough and bound fractions were collected as non-gp120-specific Ab3 and gp120-specific Ab3, respectively. These antibodies were then examined for gp120 binding activity. The data in Table 5 shows that all gp120-specific Ab3 from four monkeys bound to gp120, whereas non-gp120-specific Ab3 did

not. This finding indicates that only a fraction of Ab3 recognized gp120 and also suggests that the gp120-specific Ab3 from all four monkeys had a range of affinities for gp120. For example, gp120-specific Ab3 of PRO703 exhibited a 3 to 4-fold higher binding to gp120 than that of gp120-specific Ab3 from the other three monkeys.

Table 5
Binding of gp120-Specific Ab3 Antibodies to gp120

10	***		Binding of ¹²⁵ I-labeled <u>qp120_{SF2} to Ab3, cpm</u>				
		Ab3, ug		non-gp120 specific			
	Monkey	per well	Ab3	Ab3			
	PR0703	1	11,557	218			
15		0.2	2,324	137			
	PR0783	2	6,156	402			
		0.5	2,803	314			
	PRO419	1	4,300	236			
		0.5	1,477	NT			
20		0.25	702	NT			
	PRO431	1	5,738	212			
		0.5	2,161	NT			
		0.25	1,100	NT			

For the binding assays, RIA plates were coated with different amounts of purified anti-gp120 Ab3 for 18 hr with various amounts of antibody in PBS. The plates were blocked with 1% BSA in PBS, and 1 x 10⁵ cpm of ¹²⁵I-labeled gp120_{SF2} was added to each well. After 3 hr, the plates were washed and the radioactivities were measured. NT, not tested.

33

Neutralizing Activity of gp120-Specific Ab3

xamined the HIV-1 virus Finally, Applicant neutralizing activities of gp120-specific Ab3 antibodies. Table 6 summarizes the results. The gp120-specific Ab3 5 from PRO703 neutralized MN, IIIB and RF strains of HIV-1. The gp120-specific Ab3 from PRO431 strongly neutralized the RF strain, and weakly neutralized the MN and IIIB strains. The gp120-specific Ab3 from PRO419 and PRO783 neutralized RF and MN, respectively. The lower or limited 10 neutralizing activity of the gp120 specific Ab3s from PRO419, PRO431 and PRO783 may be due to lower affinity against certain gp120s. The gp120 binding activity of gp120-specific Ab3 from these monkeys was considerably lower than that of gp120-specific Ab3 from PRO703 (Table 5).

Table 6.

Neutralizing Activity of gp120-Specific Ab3 Antibodies

Isolated From 3C9 Immune Primate Sera

5			Neutralizing activity as % inhibition of HIV isolates		
	Monkey	Ab3, ug/ml	HIV	HIVIIIB	HIV _{RF}
	PR0703	30	97	82	84
	•	10	97	20	94
0	PR0783	55	82	0	ИЪ
-	* ***	18.8	8	0	ND
	PRO419	30	10	20	96
		15	8	o	98
		7.5	13	9	96
5	PRO431	30	46	67	96
-		15	17	56	98
		7.5	24	0	89

gp120-specific Ab3 antibodies produced by monkey PR0703

were isolated from a pool of sera obtained after the second, third, fifth, and sixth immunizations. gp120specific Ab3 antibodies from monkey PR0783 were isolated from a pool of sera obtained after the second, third, fifth, sixth, and seventh immunizations. Ab3 antibodies from monkeys PR0419 and PR0431 were isolated from a pool of sera obtained after the fourth and fifth immunizations. The neutralizing activity was assessed by comparing the amount of p24 in the well of the test antibody with that in the control well. All experiments were performed in triplicate. ND, not done.

Collectively, the results from binding and virus neutralization studies on antibodies isolated from the sera of 3C9-immunized non-human primates demonstrate that

35

3C9 can stimulate clones of primary B cells to produce br adly neutralizing HIV-1 antibodies.

Claims


- 1. A method of selecting monoclonal anti-idiotype antibody reagents useful in vaccine formulations for the treatment or prevention of HIV infection, comprising the 5 steps of:
 - (a) preparing human polyclonal anti-gp120 antibodies
 (Abls);
- (b) generating a first plurality of monoclonal antiidiotype antibodies (Ab2s) that immunologically binds with
 said human polyclonal anti-gp120 antibodies;
- (c) selecting a second plurality of monoclonal antiidiotype antibodies from said first plurality of antibodies, wherein said second plurality of antibodies are characterized by their ability to immunologically bind with anti-gp120 antibodies that neutralize in vitro multiple strains of HIV; and
- (d) selecting a third plurality of monoclonal antiidiotype antibodies from said second plurality antibodies, wherein said third plurality antibodies are characterized 20 by their ability to generate an anti-anti-idiotype antibody (Ab3) response in a primate host, wherein said antianti-idiotype antibody immunologically binds with gp120 and neutralizes in vitro multiple strains of HIV.
- The method of claim 1, wherein said human
 polyclonal anti-gp120 antibodies are pooled human sera
 from HIV positive, asymptomatic individuals.
- 3. The method of claim 2, wherein the said pooled human sera is purified for anti-gp120 antibodies by passing said pooled human sera over a column with gp120 bound to a stationary phase of said column and selecting, by elution, the bound fraction of anti-gp120 antibodies.
 - 4. The method of claim 1, wherein said first plurality of monoclonal anti-idiotype antibodies are murine.

- 5. The method f claim 1, wherein said anti-gp120 antibodies of step (c) are human.
- 6. The method of claim 1, wherein said anti-gp120 antibodies of step (c) are human antibodies selected from 5 said human polyclonal antibodies of step (a).
- 7. The method of claim 6, wherein said anti-gp120 antibodies of step (c) are selected by passing said human polyclonal antibodies of step (a) over a column with an anti-idiotype antibody from said first plurality of antibodies bound to a stationary phase of said column, eluting bound fraction of anti-gp120 antibodies and selecting from said bound fraction those anti-gp120 antibodies that are capable of neutralizing in vitro multiple strains of HIV.
- 8. A monoclonal anti-idiotype antibody selected according to the method of claim 1.
 - 9. A monoclonal anti-idiotype antibody that is capable of eliciting a group neutralizing anti-HIV antibody response in humans.
- 20 10. A hybridoma cell line capable of producing a monoclonal anti-idiotype antibody that is capable of eliciting a group neutralizing anti-HIV antibody response in humans.
- 11. A hybridoma cell line deposited with the ATCC, 25 Rockville, MD, under accession No. HB 10701.
 - 12. A monoclonal antibody produced by the cell line of claim 11.
 - 13. A monoclonal antibody which elicits a substantially equivalent immune response in humans as a


monoclonal antibody produc d by the hybridoma cell line of claim 10 or 11.

- 14. A monoclonal antibody which reacts with and competes for binding to the same anti-gp120 antibody as a monoclonal antibody produced by the hybridoma cell line of claim 10 or 11.
- 15. A vaccine formulation comprising, in admixture with a pharmaceutically acceptable carrier, a monoclonal anti-idiotype antibody that is capable of eliciting a group neutralizing anti-HIV antibody response in humans.
 - 16. A vaccine formulation comprising, in admixture with a pharmaceutically acceptable carrier, an antibody of claim 8 in a dose sufficient to induce an active immune response in humans.
- 17. A vaccine formulation comprising, in admixture with a pharmaceutically acceptable carrier, an antibody of claim 9 in a dose sufficient to induce an active immune response in humans.
- 18. A vaccine formulation comprising, in admixture 20 with a pharmaceutically acceptable carrier, an antibody of claim 12 in a dose sufficient to induce an active immune response in humans.
- 19. A vaccine formulation comprising, in admixture with a pharmaceutically acceptable carrier, an antibody of claim 13 in a dose sufficient to induce an active immune response in humans.
- 20. A vaccine formulation comprising, in admixture with a pharmaceutically acc ptable carrier, an antibody of claim 14 in a dose sufficient to induce an active immune 30 response in humans.

A panel of monoclonal anti-idiotype antibodies

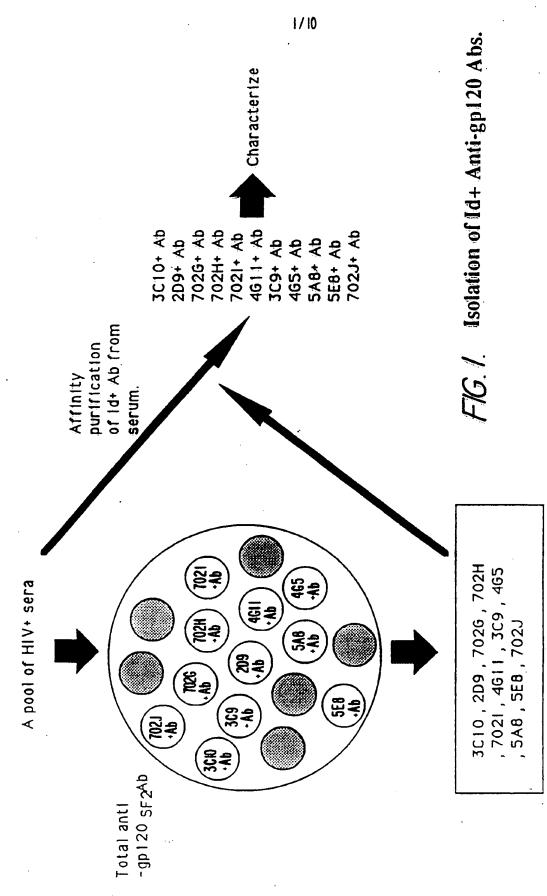

SUBSTITUTE SHEET

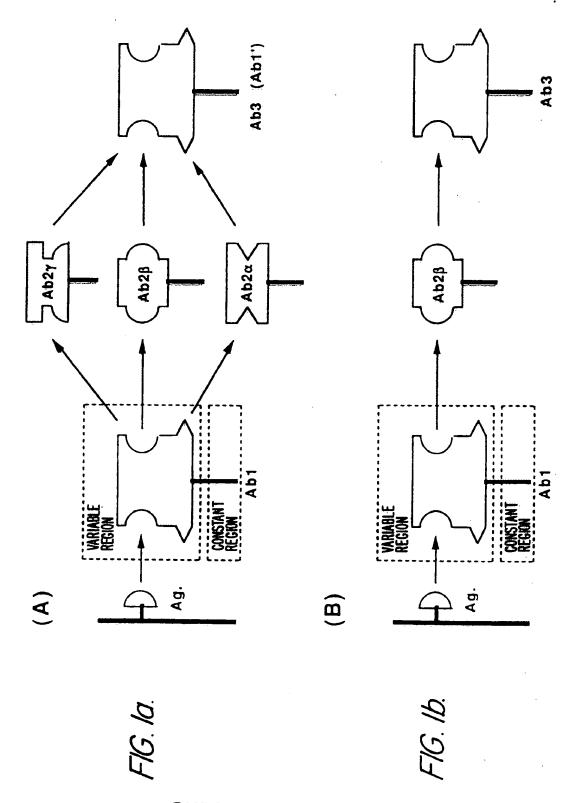


BUBSTITUTE SHEET

3/10


F/G. 2.


00001


.001

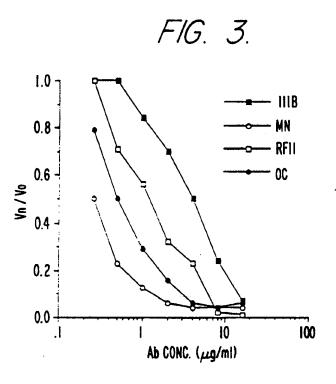


SUBSTITUTE SHEET

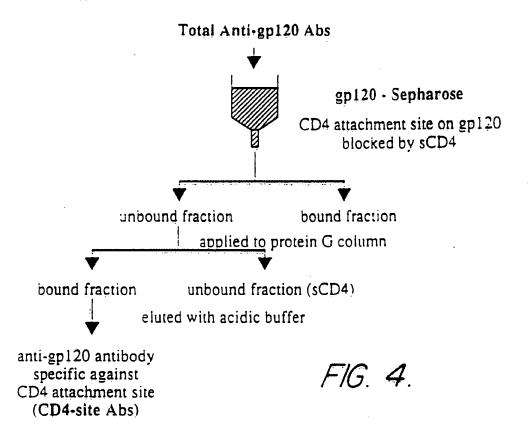

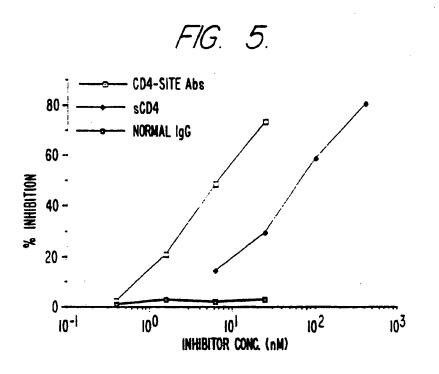

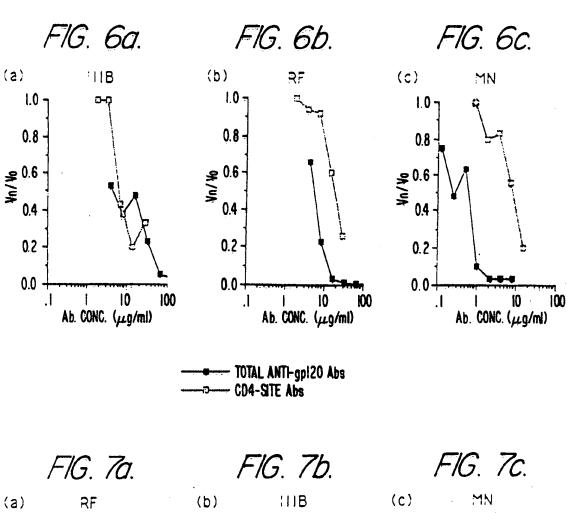

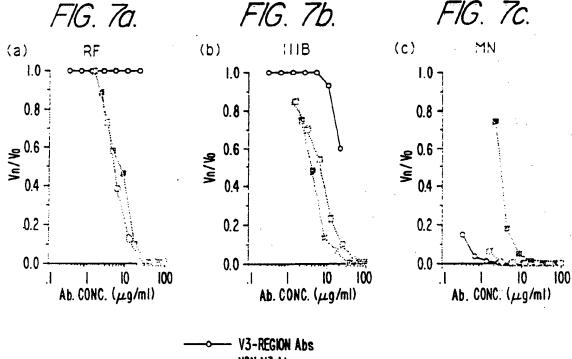

NON-V3 Abs
TOTAL ANTI-gp 120 Abs

HIVRF Virus Neutralization In HIVIIIB Ab1s HIVMN Corresponding Anti-id Mabs +++ +++ 02 Characterization of the Against Different FIG. 8. Binding to gp120rms Characteristics of Abls Binding to gp120SF2 1A12 id 2A12 id 4G11 id 4E10 id 2D9 id 3A6 id 1C9 id 3C9 id 2F2 id 5A8 id 4G5 id SE8 id Abls Not done Ab2s 1 A 1 2 3C10 2 A 1 2 4E10 ND: 4611 316 2 D 9 1C9 3C9 2 F 2 465 5 A 8

parenthesis

Characterization of the Corresponding Ab3s Against Different Ab2 in Monkeys

		Binding to		Virus Neutralizations in	ns in
A b 2 s	A b 3 s	gp 1 2 0	HIVMN	HIVITIB	HIVRF
3C9	3C9 Ab3	++ (2/2)*	+	+	+
4 G 1 3	4G11 Ab3	. (0/2)	•	•	•
4 (3.5	4G5 Ab3	+ (2/2)		,	•
5 A 8	5A8 Ab3	. (0/2)		•	•


*Number of positive responses from number of immunized monkeys in


Neutralizing Activity of Ab3 Isolated From Monkey Sera

•	7 4 V		as % Inhibition in	in
5		1	gran and	
9 Ab3*	30 µg/ml	9.1%	82%	84%
	10 µg/ml	97%	20%	94%


*Purification was done on Ab2 and gp120 affinity column

FIG 10

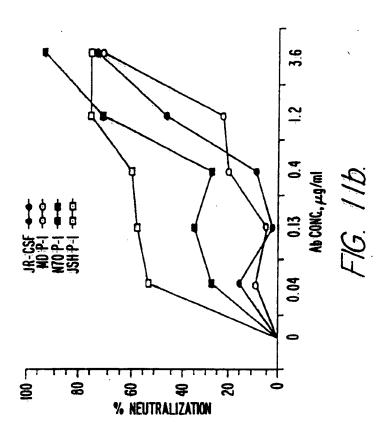

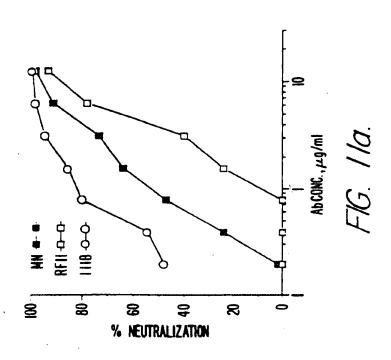

SUBSTITUTE SHEET

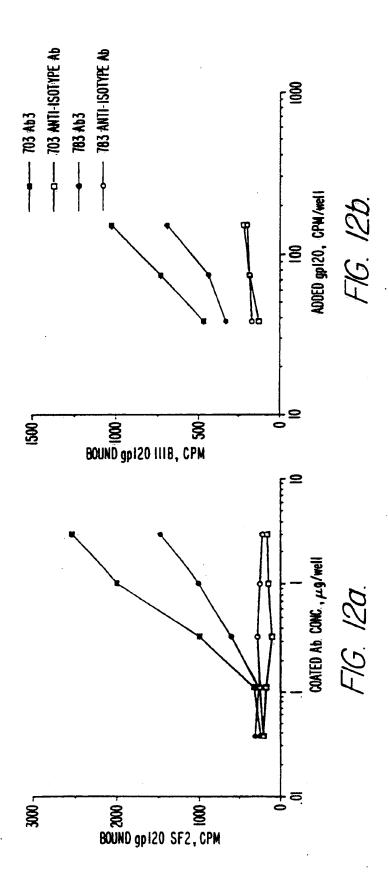

SUBSTITUTE SHEET

1			t (If several classification				
	international Patent 5 GO1N33/5		GO1N33/68;		рс (39/21;	C12	P21/08
II. FIELDS SI	EARCHED		The second secon				
			Мівівня Восиц	estation Searched	,		
Classification	System		and the same of the state of the same of t	Classification Syr	a bols		
Int.Cl.	5	GOIN;	A61K				
		D to the	ocumentation Searched other Extent that such Documents	r than Minimum D are Included in th	ocumentation • Fields Searche	48	
III. DOCUME	NTS CONSIDERE						
Category o	Citation of De	equent, 11 m	ith indication, where approp	riste, of the relevan	ot passages 12		Relevant to Claim No.13
X	vol. 36 CANADA pages 8 R.C. KE	, no. 1: 11 - 810 NNEDY E	AL OF MICROBIOL 1, 1 November 1 5; 1 AL.: 'Concept epatitis b viru	990, MONTR s of idiot	ype-based	1	1,4-14
Y	immunodi see pagi column	eficience 814, c 2, line	cy virus. column 2, line 17; figures 1, (G.B. THORTON	49 - page 3	815,		1-20
		umn 6,	line 16 - colum	n 8, line -/	5		
		10		opo leter dece	ment published a	fter the internal	ional filing date
"A" docum consid "E" earlier filing	date	peral state of plar relevance lished on or a	ter the international	or priority cited to u invention "X" document cannot be	of particular role considered novel in the print of particular role considered novel in inventive step	conflict with the neighe or theory evence; the cisis	application but underlying the sed invention
which citatio "O" docum other	means	the publications on all disclosur	no date of another	'Y' document cannot be document ments, su in the art	of particular rele considered to in is combined with the combination b	volve as inventi- heing obvious to seing obvious to	re step when the her such docu- a person skilled
	than the priority dat	e cialmed		. 4. 40cm	member of the s	erne berein rem	.,
IV. CERTIFIC		the Internal	nal Sannch	Date of M	alling of this inte	emational Search	h Report
mate of the Ac	tual Completion of	JUNE 19			06. 92		Î
International S	earching Authority			1 -	of Authorized Of	0	10 300
	EUROPE	AN PATEN	T OFFICE	VA	N BOHEMEN	N C.G.	19 ran Belina

	NTS CONSIDERED TO BE RELEVANT (CONTINUED FROM THE SECOND SHEET)	
	ONSIDERED TO BE RELEVANT (CONTINUED FROM THE SECOND SHEET) Citation of Document, with indication, where appropriate, of the relevant passages	Relevant to Claim No.
ategory °	Citation of Document, with indication, where appropriate, at	
		1 00
	AIDS RESEARCH AND HUMAN RETROVIRUSES vol. 6, no. 1, 1 January 1990, NEW YORK NY USA	1-20
	page 42;	
	antibodies to a neutralizing antibody (bat123) recognizing an immunodominant domain of hiv-1	
	gp120.' see the whole document	
	त शंच	
	•	
	·	

ANNEX TO THE INTERNATIONAL SEARCH REPORT ON INTERNATIONAL PATENT APPLICATION NO. US 9201749 58376 SA

This annex lists the patrat family members relating to the patent documents cited in the above-mentioned international search report.

The members are as contained in the European Patrat Office EDP file on
The European Patent Office is in an way liable for those sarticulars which are merely given for the purpose of information. 23/06/92

Patent document cited in search report	Publication date	Patost family monitor(s)	Publication date
US-A-4908203	13-03-90	None	
•			
·			
			•
•			
more details about this annex :			